Certification



State Food and Drug Administration



NRTL



CE (MDD)





ISO13485 / ISO9001

Sleeves

Picture	Sleeve Type	Order NO.	Chamber
	Arm Sleeve	220041	4
-		220081	8
	Leg Sleeve	220042	4
		220082	8
	Waist Sleeve	220083	4
		220083L	
homes of	Full Body Sleeve	220016FB	8
-1	Calf Sleeve (DVT Sleeve)	220052	2
	Foot Sleeve (DVT Sleeve)	220053	1



Guangzhou Longest Science & Technology Co., Ltd.

Address: #96 Chuangqiang Road, Ningxi Street, Zengcheng District, Guangzhou, China 511399

Tel: 86-20-66353999 Fax: 86-20-66353920

E-mail: export@longest.cn **Website:** www.gzlongest.com








Health Life Longest Care

Longest




• Mechanism and functions

The inflation and deflation of the sleeves enhances the backflow of venous blood and lymph, which can improve the blood circulation of diseased areas. This will lead a series of beneficial effect such as improving the function of peripheral blood vessels, removing the edema, accelerating the healing process, and preventing Muscular Dystrophy.

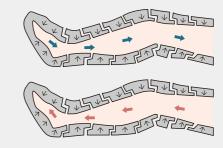

Features

- Overlapping air chambers generate a gliding pressure wavewith high efficiency
- ① User-friendly interface, just a few steps to finish your setting
- O Gentle inflation and deflation
- Real-time pressure, time, and inflation part monitoring (visible in the interface)
- Pressure unit switchable between mmHg and kPa

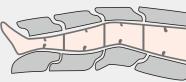

Overlapping Air Chamber

Four chambers modes





-our chambers modes

• Specifications

Channel	2
Sleeve chambers	1/2/4/8
Compression modes	18 (Including 3 DVT modes + 6 four chambers modes+ 9 eight chambers modes)
Time range	1 - 99 m1ns
Pressure range	37 - 187 mmHg
Interface	5.7" touch screen
Dimension	393 (L) x 366 (W) x 194 (H) mm
Weight	10kg

Application

- Prophylaxis of Deep vein thrombosis (DVT) (Post-stroke deep vein thrombosis, Post-traumatic deep vein thrombosis)
- ② Prophylaxis and treatment of edema (Chronic venous edema, Lymphedema, Post-mastectomy lymphedema



